



Manual despiker for spin data. (Work in progress, 2025.11.07)

```
data = dopey.load("file.xy")
dopey.plot(data)
```

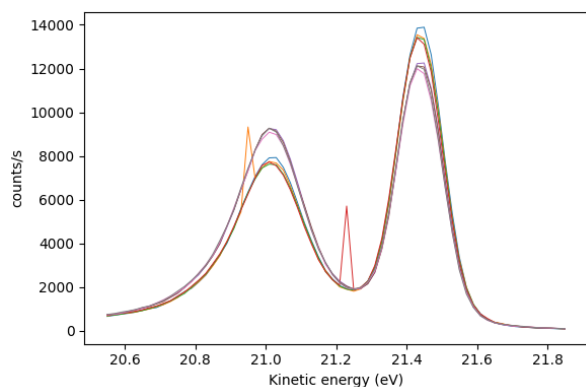


Pass the data into the manual despiker method. It currently works for spin EDC and spin MDC. Select a curve with a spike with the *Curve slider* and the point with the *Point slider*:


```
dopey.despikeSpinManual(data)
```

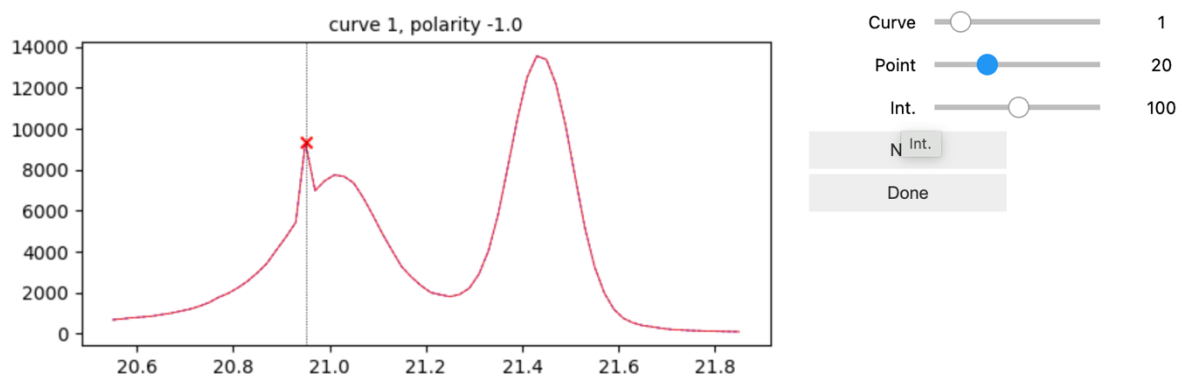


Adjust the intensity in that point with the *Int. slider*:

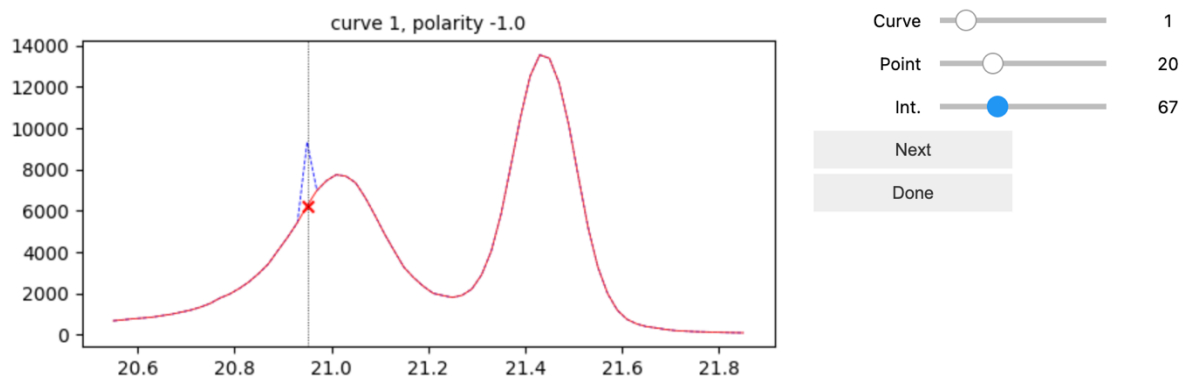

Click the Next button when ready and, if necessary, select another curve with the *Curve slider* and repeat:

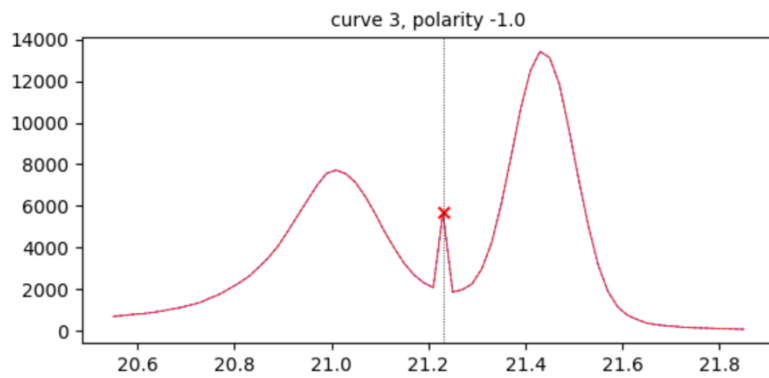


Curve 3
 Point 34
 Int. 100
 Next
 Done

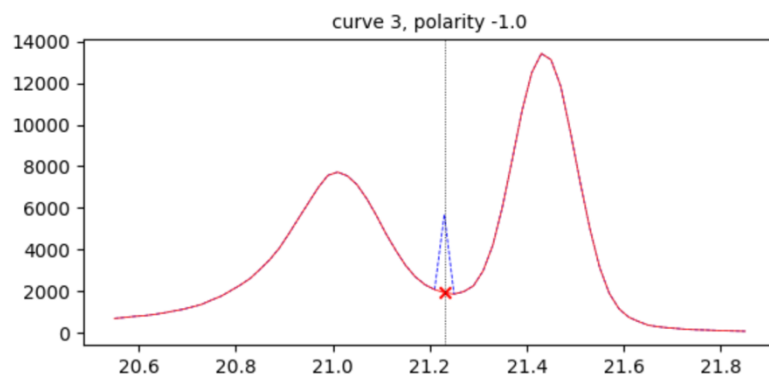

Curve 3
 Point 34
 Int. 34
 Next
 Done

Click the *Done* button when ready and plot the data.

`dopey.plot(data)`

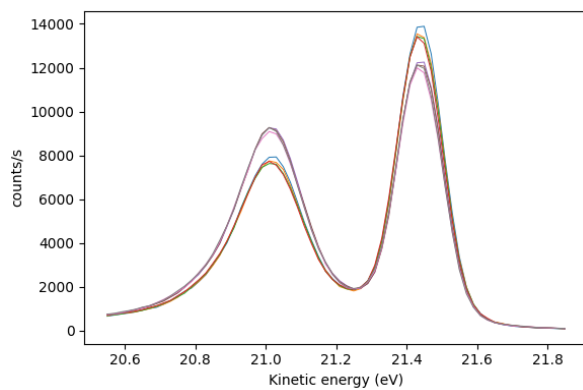

Save the data to file so that this does not have to be repeated. Use the `pickleSave()` method:

`dopey.pickleSave(data, "data.pickle")`

Saved to file `data.pickle`

The data can be loaded with `pickleLoad()`:

`data = dopey.pickleLoad("data.pickle")`

Loaded file `data.pickle`
 Data type `spin_edc`